EN

DE

CZ

SK

PL

HU

SI

RO

UA

RU

SONIC TS400

USER MANUAL  
TESLA TOOTHBRUSH

## PRODUCT OVERVIEW

ULTRASONIC WHITENING TOOTHBRUSH  
5 BRUSHING MODES  
3 LEVELS OF INTENSITY

## IMPORTANT INSTRUCTIONS

### General notices

- This device is not intended for use by persons (including children) with reduced physical sensory or mental capabilities.
- Children between the age of 3-10 years must be given supervision or instruction concerning use of the device by a person responsible for their safety. It is not recommended for use of children under the age of 3 years.
- This device is only for personal use. It not intended for use in dental clinics by multiple patients.
- This device is for cleaning teeth, gums and tongue only.
- Charge indicator may not work instantly if battery fully discharged. It will switch back on in 6-10 minutes.
- It is normal to feel slight heat that comes through the handle during charging process. If it turns extremely hot, turn the charger off immediately.
- Store this device away from direct sunlight at temperature at 7 °C ~ 40 °C- Do not expose surface of this device to certain substance (gel, oil, soap etc.) to extend lifespan and to avoid corrosion.
- Toothbrush heads used daily and for a long time leads to a large growth of microbes and bacteria between the bristles, which is not good for cleaning the mouth. it is recommended to replace the toothbrush heads every 3 months.

### Warning

- Do not place the charger above or near the water source. For example a bathtub, a pool etc.
- Do not change batteries inside the device.
- It is normal for first-time users to experience slight bleeding and discomfort for their initial week of use. If the bleeding or discomfort continues to occur after 1 week of use contact your dental professional.
- If you had an oral or gum surgery in the last 2 months or implanted dental prosthesis, consult your dental professional prior to use.


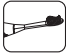
## CHARGING INSTRUCTIONS

- If the charging indicator is flashing, it needs to be charged.
- Plug the charging cable to any USB outlet then set the device onto the charging base to start charging. The indicator lights up (up to 4 charged levels).
- When fully charged, the device can be used for up to 180 days.


## MODES

- Press the power button, the device will operate in the mode it was last in before turning off.
- If you want to change the mode press the power button, you will have 5 s to select the desired mode. If you fail to select the desired mode within 5 seconds, you must turn the device off and on again and select the desired mode within 5 seconds.
- To turn the device off, press and hold the power button for 2 s.

## GETTING STARTED

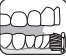	<p>Insert the brush head to the drive shaft while a certain gap must be left between for the normal vibration.</p>		<p>Place the bristles against the teeth at 45o towards the gum line.</p>
	<p>Wet the bristles and apply a certain amountt of toothpaste.</p>		<p>Move the bristles slowly in circular motion. Do not scrub the teeth, but apply light pressure to let the device do the brushing for you.</p>
	<p>Put it into your mouth before switching the device on to avoid accidental toothpaste splash. Then press the power button again to select the mode.</p>		<p>After brushing your teeth, gargle your mouth.</p>

## BRUSHING TIME AND INTENSITY CONTROL

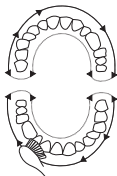The timer features 30 second intervals for the full dentist recommended 2 minutes of brushing. There is a short pause after each interval reminding you to switch to a different area.

5 cleaning modes.  
3 levels of intensity: weak - medium - strong.  
To select vibration speed press intensity button multiple times to set desirable speed.

## FUNCTION

2 min timer

180 days standby time

IPX7 waterproof

Wireless charging

5 cleaning modes  
3 levels of intensity = 15 variation to select

## MAINTENANCE

### Clean the brush head

After each cleaning, rinse the brush to remove residue and allow it to air dry.

Do not rub the bristles to avoid deforming them.

### Clean the handle

Rinse it after each brushing to remove any residual debris.

Do not immerse the device in water for a long time.

Do not clean the toothbrush with water at a temperature higher than 40 degrees or with a corrosive cleaning agent.

### Replace the toothbrush head

It is recommended to replace your toothbrush head every 3 months for the best brushing results.

Also replace the toothbrush head when the bristles are bent outwards.

## TROUBLESHOOTING

Problem	Diagnosis	Solution
Weakened vibrations.	<ol style="list-style-type: none"><li>1. The bristle is pressed too hard</li><li>2. The brush head is fixed too tightly</li></ol>	<ol style="list-style-type: none"><li>1. Do not press the bristles with too much force.</li><li>2. Reconnect the brush head to the drive shaft and make sure it is fixed with some gap between the brushes.</li></ol>
The toothbrush can only work for a short period of time.	<ol style="list-style-type: none"><li>1. The toothbrush has not been used for more than 6 months.</li><li>2. Incorrect charging</li><li>3. Insufficient charging time</li></ol>	<ol style="list-style-type: none"><li>1. Make sure the device is placed on the charging seat in an upright position and the charging pin is inserted.</li><li>2. Charge properly.</li></ol>
The toothbrush does not work at all.	<ol style="list-style-type: none"><li>1. Failure of internal mechanics</li><li>2. Battery almost reached its lifetime</li></ol>	Send for complaint.

## SPECIFICATION

Vibration Frequency:	25500-48000 strokes/min
Input parameters:	DC 5 V, 1 A
Out parameters:	DC 3.7 V, 300 - 500 mA
Charging time:	24 hours
Battery:	Lithium Battery
Working time:	About 180 days (4 mins/day)

### Package content:

Tesla Toothbrush Sonic TS400  
2x Toothbrush heads  
2x Hygienic caps for toothbrush heads  
1x Wireless charging base  
1x Travel case  
Manual

## INFORMATION ABOUT DISPOSAL AND RECYCLING

This product is marked with the symbol for separate collection. The product must be disposed of in accordance with the regulations for disposal of electrical and electronic equipment (Directive 2012/19/EU on waste electrical and electronic equipment). Disposal together with regular municipal waste is prohibited. Dispose of all electrical and electronic products in accordance with all local and European regulations at the designated collection points which hold the appropriate authorisation and certification in line with the local and legislative regulations. Correct disposal and recycling help to minimise impacts on the environment and human health. Further information regarding disposal can be obtained from the vendor, authorised service centre or local authorities.

## EU DECLARATION OF CONFORMITY

Hereby, Tesla Solar s.r.o. declares that the radio equipment type TSL-PC-TS400 is in compliance with 2014/53/EU directives. The full text of the EU declaration of conformity is available at the following internet address: [tsl.sh/doc](http://tsl.sh/doc)

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